

Training Manual



Debussy
Thermodermie®

*B*eauty from Within

Contents

Part 1 - Understanding Thermodermie®

p.3

- 1.1 : Thermodermie in action ? p.4/5
- 1.2 : Results on the skin p.6
- 1.3 : Debussy® method p.7

Part 5 - Debussy® Facial

p.28

- 5.1 : Description of the treatment p.29
- 5.2 : The Protocol p.30
- 5.3 : The Steps p.31
- 5.4 : The Movements p.32

Part 2 - Beauty face® Presentation

p.8

- 2.1 : The machine p.9
- 2.2 : The Control Panel p.10
- 2.3 : Maintenance p.11
- 2.4 : Trouble shooting p.12/13
- 2.5 : Safety Guidelines p.14

Part 6 - Debussy® Vitality & Rejuvenation

p.33

- 6.1 : Description of the treatments p.34
- 6.2 : The Protocol & Steps (Session 1) p.35
- 6.3 : The Protocol & Steps (Sessions 2,3,4 & 5) p.36
- 6.4 : The Movements p.37

Part 3 - Using Debussy® Accessories

p.15

- 3.1 : Presentation of the Accessories p.16
- 3.2 : The Facial Tube (Suction & Infrared) p.17/18
- 3.3 : The Ultrasound Probe (UP) p.19/20
- 3.4 : The Beauty Cleaner® p.21/22

Part 7 - In the Booth

p.38

- 7.1 : Setting up the booth p.39
- 7.2 : Keeping your clients informed p.40
- 7.3 : Your treatment worksheets p.41/42/43/44

Part 4 - Beauty Cleaner®

p.23

- 4.1 : Description of the treatment p.24
- 4.2 : The Protocol p.25
- 4.3 : The Steps p.26
- 4.4 : The Movements p.27

Part 8 - Cosmetic Line

p.45/46

Part I

Understanding Thermodermie®

I.1 : Thermodermie in action ?

I.2 : Results on the skin

I.3 : DEBUSSY® method




I.1 : Thermodermie® in action ?



Beauty Face®




Innovative equipment which combines
4 cutting edge technologies:

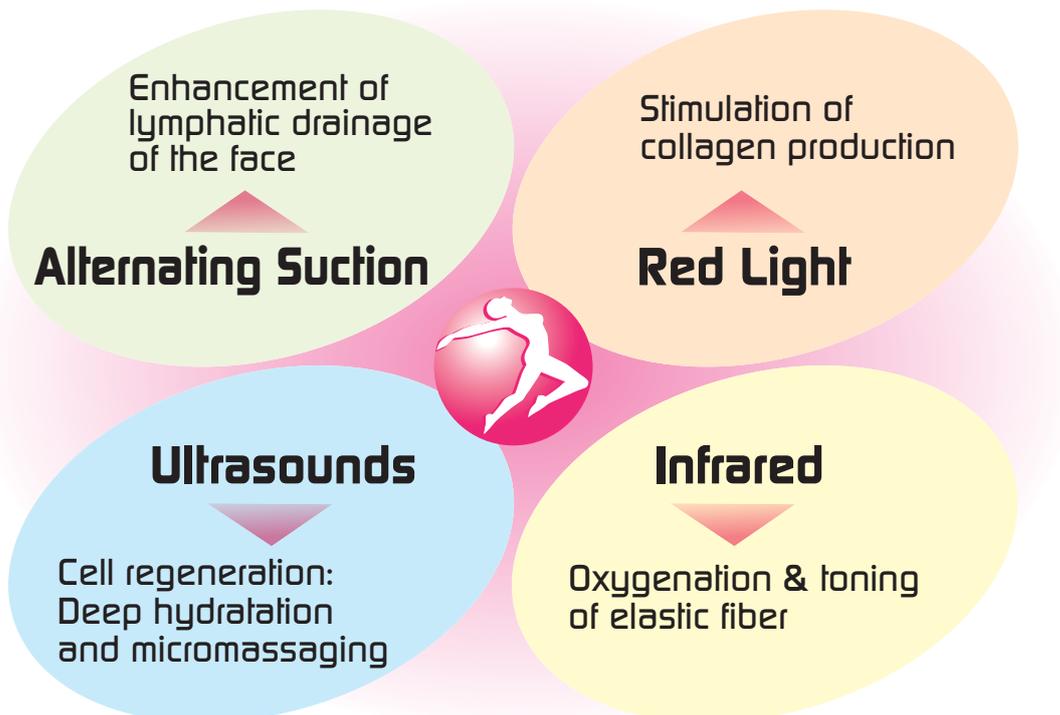




I.1 : Thermodermie® in action ?

Restoring the face's own resources,
slowing down the natural ageing process.

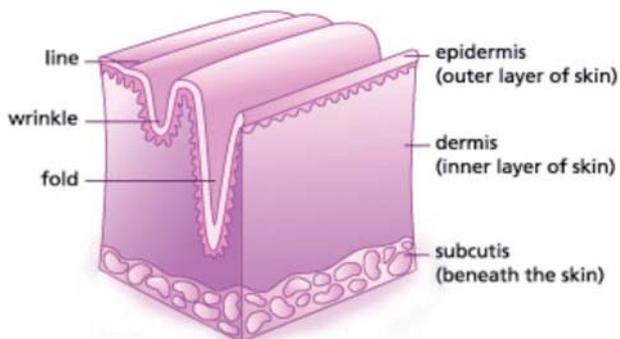


Simultaneously acts on
all 3 skin layers.

Ageing collagen

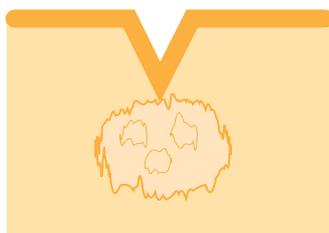

Loosening of the skin, dull
complexion, visible wrinkles.

Thermodermie®
treated collagen

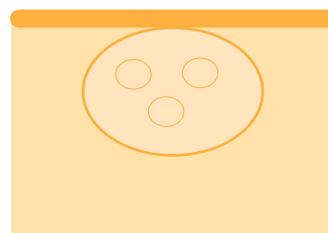

After treatment, the skin
is tighter and brighter.

1.2 : Results on skin

Fresher, Healthier, Younger skin.

Deep cleansed skin

- ▶ Detoxified and oxygenated tissues
- ▶ Enhanced lymphatic system

Increased collagen production

- ▶ Firmer skin
- ▶ Remodelled contour of the face

Deep hydrated skin

- ▶ Toned elastic fibres
- ▶ Tightened pores



I.3 : DEBUSSY® Method

Non-Invasive

- ▶ No injection, no operation
- ▶ Relaxing massage of the dermis

Fast

- ▶ Your client can choose:
1, 3 or 5 sessions.
- ▶ Full treatment can be done
in only 5 days (7h30).

Effective

- ▶ Visible results after first session
- ▶ Long lasting fresh look

- ▶ Clinically and dermatologically tested.
- ▶ Debussy 's advice : drink mineral water
(minimum 1,5 liters) during the treatment.



Part 2

Beauty face® Presentation

2.1 : The machine

2.2 : The Control Panel

2.3 : Maintenance

2.4 : Trouble shooting

2.5 : Safety Guidelines



Beauty Face®

2.1 : The Machine



Control panel



Tube and nozzle
for facial massage:
suction and infrared

Casing:
includes motor

2.2 : The Control Panel



Installation block
for suction tube

Plug for infrared

Power switch : ON/OFF


■ Press the red indicator: ON/UP
(ultrasound probe)

■ Press the orange indicator:
alternating suction (sensitive skin)

Usage meter
in hours

Plug for ultrasound
probe (UP)

Power knob:
to increase suction
(4 = maximum)

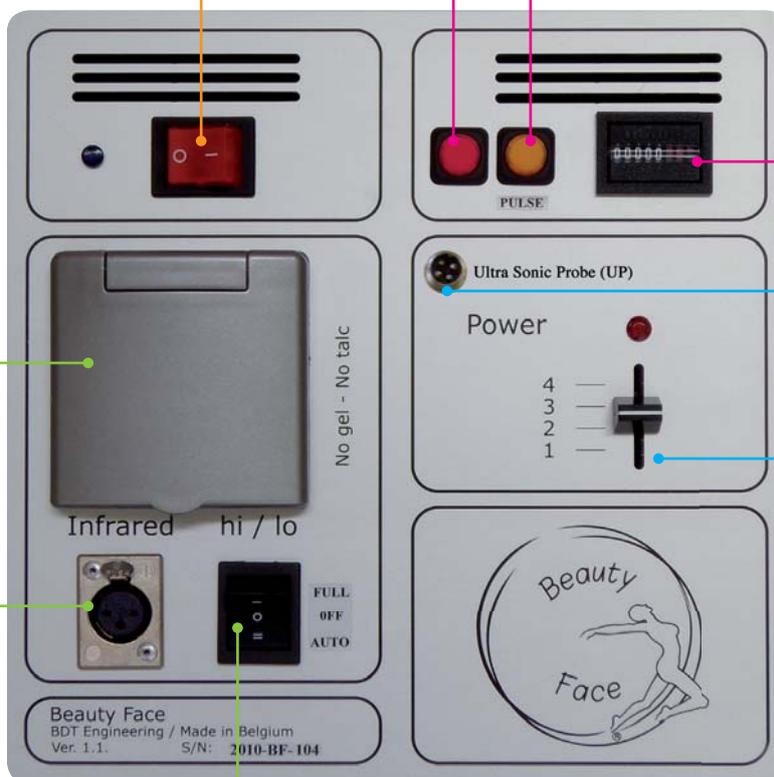

Infrared adjustment knob:
Full (-) : maximum
Off (0) : inactive
Auto (=) : automatic/variable

2.3 : Maintenance



THE BEAUTY FACE®

Must be cleaned twice a day. In the morning before the first treatment and at the end of the day after the last treatment. Turn off the machine before cleaning and use a dry cloth. Never use: thinners, detergents, acidic solutions and flammable liquids.



FACIAL NOZZLE
AND ULTRASOUND
PROBE (UP)

Must absolutely be cleaned after each treatment with an antiseptic wipe and dry cloth.



THE INFRARED
BULBS

Dust can build up behind the bulbs. You must unscrew the bulbs and clean regularly.
The surface of the bulb must be wiped each day with a dry cloth.
Never wipe a hot bulb.
Never use alcohol.
Never touch the lamp with bare fingers.



THE FILTER

Must be cleaned once a day. Remove from the suction block and rinse with antibacterial soap. Make sure to dry the filter before placing it back in the suction block.

Important : if you have a sterilization machine, it is strongly urged to use it, once a day, for all the Beauty Face® detachable accessories.

2.4 : Trouble shooting



TO CHECK IF THE ULTRASOUND PROBE FUNCTIONS PROPERLY:

Put a little water on the steel surface of the ultrasound probe. If the water evaporates and there are small bubbles, then the ultrasound probe is working properly. Otherwise, the ultrasound probe is defective or broken. You must then replace it and order a new ultrasound probe through your local Debussy supplier (Bali).


WHAT IF THE MOTOR STOPS?

- ▶ Due to over heating, the machine will stop automatically. Wait 15 minutes and the machine will restart normally.
- ▶ Check that the power cord is correctly plugged in and that the wall socket is working properly. If there is a problem, contact an electrician.
- ▶ A fuse may have melted. Replace it (see next page)

WHAT IF...?

- ▶ A bulb does not work: change it by unscrewing the facial nozzle. Infrared bulbs have a lifespan of 4000 to 5000 hours.
- ▶ Suction is weak: check that dust has not built up behind the bulbs (clean) or in the filter (clean).

2.4 : Trouble shooting

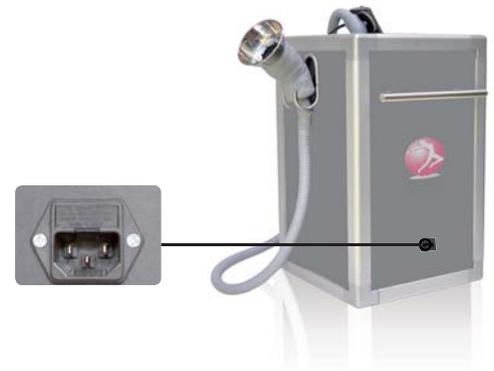

The fuses are located behind the Beauty Face®, just above the electric connection of the machine. It is a black little black box. (cf. Sketch to change the fuse)




Pull out the drawer with a small screwdriver.



Remove the used fuse then the new one.




Put the next fuse instead of the used one.

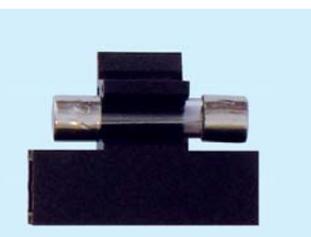

Eventually, put a new fuse in the reserve.

2.5 : Safety Guidelines

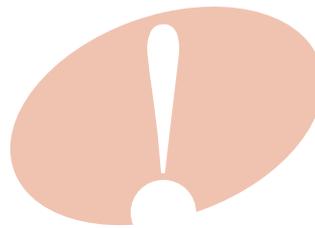

1

You must neither obstruct the ventilator's air passage nor open and take apart the machine. If there is a problem, contact your manager.

2

The ultrasound probe is very expensive and is sensitive to shocks. It is therefore necessary to use it with caution and to never drop it.

3

Follow closely the cleaning instructions of the Beauty Face®. You must not use: thinners, detergents, acidic solutions and flammable liquids.

4

The infrared light is hot. While giving treatment, inform and question the client. Do not hesitate to personally check the heat level on the skin of your hand.

5

The Beauty Face® must be turned off : after each treatment, before being cleaned and before any form of maintenance or repair. At the end of the day, unplug the Beauty Face® and put away all the accessories.

Part 3

Using DEBUSSY® Accessories

3.1 : Presentation of the Accessories

3.2 : The Facial Tube (Suction & Infrared)

3.3 : The Ultrasound Probe (UP)

3.4 : The Beauty Cleaner®



3.1 :

Presentation of the Accessories



Ultrasound Probe (UP)

Plug into Beauty Face®.

Purpose :

- ▶ Tissue hydration
- ▶ Cell massage

Tube and facial nozzle

Plug into Beauty Face®.

Purpose :

- ▶ Suction : lymphatic drainage
- ▶ Infrared: collagen stimulation



Facial
nozzle

Infrared
50W

Beauty Cleaner® (Ultrasound)

Purpose :

- ▶ Exfoliation (deep cleaning and soft peeling)
- ▶ Balancing of the skin's PH level



3.2 :

The Facial Tube

(Suction & Infrared)



INSTALLATION OF THE FACIAL TUBE (suction and infrared):

- 1 Open the gray hatch to install the end of the facial tube.
- 2 Then attach the 3 prong plug (infrared plug) to the black plug.
- 3 To unplug the infrared plug press on the « PUSH » button and pull.


SETTINGS TO FOLLOW:

- Normal skin = HI/LO : FULL + POWER : 1 + continuous suction (orange button « Pulse » not pressed down).
- Sensitive skin = HI/LO : FULL + POWER : 1 + pulse suction (orange button « Pulse » pressed down).



INFRARED SETTINGS HI/LO:

- FULL (-) : maximum
- OFF (0) : inactive
- AUTO (=) : auto/variable

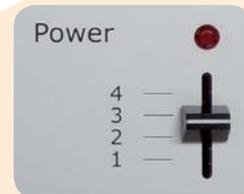

USING THE «POWER» SWITCH:

The POWER switch is used to increase suction power (massage). Level 4 is maximum power.



USING THE «PULSE» SWITCH:


The pulse switch is used to alternate between suction and infrared.

- Orange button not pressed down = continuous suction and infrared.
- Orange button pressed down = pulse mode. TO USE ON THIN OR VERY SENSITIVE SKIN.

3.2 :

The Facial Tube

(Suction & Infrared)

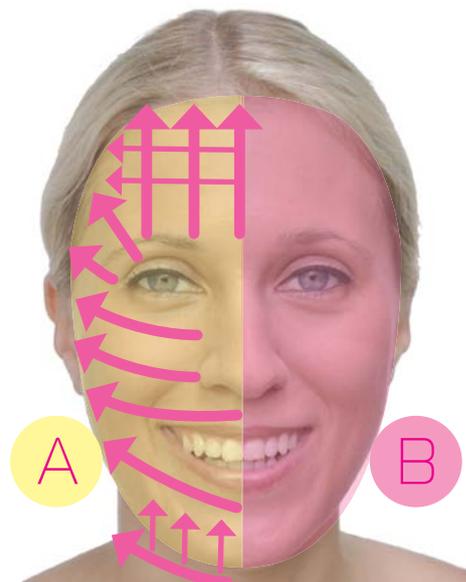

The massage is a two-step process:
A left side and B right side

USING THE FACIAL TUBE (Suction + Infrared)

- ▶ The facial massage consists of a series of even lines covering the surface of the face (see part 5, treatment protocols). Starting from the neck and progressing steadily to the forehead.
- ▶ Each line starts towards the interior of the face and continues towards the exterior (from chin to lower ear, from nose to the hairline).
- ▶ In order to optimize the tightening effect, each line rises slightly towards the end of the movement.
- ▶ The facial nozzle must systematically be pressed onto the skin and not allow any air through.
- ▶ One hand uses the tube for the massage while the other hand precedes the movement and stretches the skin in front of the passing tube.
- ▶ When wrinkles are visible, you must stretch the wrinkles with your fingers in front of the passing facial tube.
- ▶ The infrared light warms up. During the treatment, you must both monitor the heat level on the face with your fingers and ask the client about their sensations.

3.3 :

The Ultrasound Probe (UP)



INSTALLING THE “UP”:

Plug the Ultrasound Probe into the three-pronged black socket (located above «Power»).



SETTING GUIDELINES:

Press the POWER button on the probe to turn it on.

- ▶ Low : Normal skin.
- ▶ Pulse : Sensitive skin.


After 10 minutes of use, a beep will signal that the probe is turning off (total time of the treatment for the whole face).

3.3 :

The Ultrasound Probe



10 min

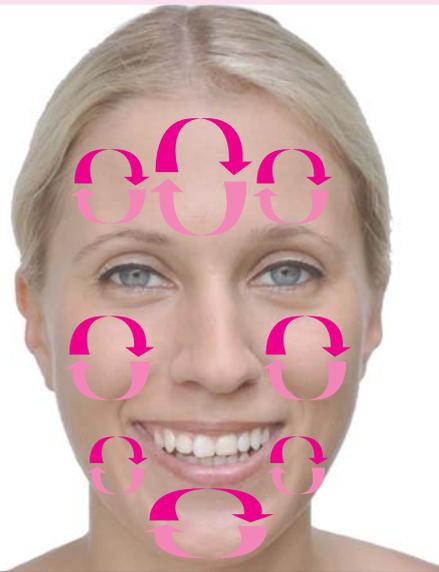

10-minute program to aid the aloe vera gel to penetrate throughout the skin of the face.

USING THE ULTRASOUND PROBE

- ▶ When using the probe, make sure to cover evenly the skin with the ALOE VERA gel.
- ▶ The steel head must be flat and in permanent contact with the skin. You should never hold the probe slanted. Any air passing through is to be avoided.
- ▶ Perform slow and circular clockwise movements over the whole face, only slight pressure is needed (it is not necessary to push down).

3.4 :

The Beauty Cleaner®



INSTRUCTIONS

- 1 Remove the plastic protective tip.
- 2 Turn on the machine by pressing the ON/OFF button (red LED turns on = working).
- 3 Press the VOL button (red LED), Select the mode Low, High, with or without W (pulse mode).
- 4 The machine turns off automatically after 10 minutes (emits 2 beeps).

POWER MANAGEMENT

- When the battery is low, the machine turns off after emitting 5 beeps.
- When charging : three blinking LED (green, orange, blue).
- When charging is complete, the three LED stop blinking and stay on.


SETTINGS

Press the VOL button (red LED) and you can select in order from the following modes:

- L (Low, green LED) = low power mode. For sensitive skin or for delicate areas (under the eyes).
- H (High, blue LED) = high power mode. For normal skin.
- W (Pulse, orange LED) = pulse mode. Can be used with modes L or H; when targeting a specific area (under the eyes, acne). This mode helps to avoid overheating of the skin.
- P (red LED) = in use.


3.4 :

The Beauty Cleaner®

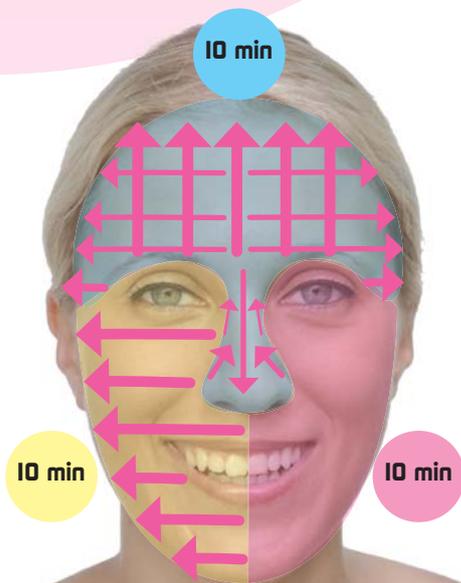

Treatment Beauty Cleaner® (alone)

or

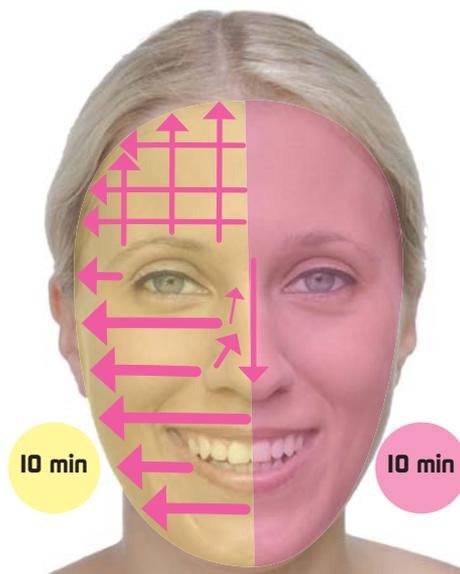

Treatment Beauty Cleaner® with Beauty Face®

USING THE BEAUTY CLEANER

The skin must **IMPERITIVELY** and **CONTINUOUSLY** be moistened (mineral water) when using the **BEAUTY CLEANER®**.

You use a moist sponge or spray before each facial movement of the **BEAUTY CLEANER®**. The skin is moist enough, when you can see a little water vapor in front of the head of the machine.

- ▶ The head of the **BEAUTY CLEANER®** must be wiped at regular intervals during the treatment to avoid deposits. Maintain the handle of the machine at a 45° angle on the skin so that the head is vertical (see photo).
- ▶ It is not necessary to push down the head onto the skin. The head should «surf» over the curves of the face.
- ▶ The movements always start from the center of the face and continue towards the exterior (see part 4 treatment protocols). Perform even lines starting from the oval of the face and progressing towards the forehead. Do not forget to work on the nose. For a good treatment, three passages minimum are necessary over the same part of the face.

Part 4

Treatment Protocols

Debussy® Beauty Cleaner®

4.1 : Description of the treatment

4.2 : The Protocol

4.3 : The Steps

4.4 : The Movements



4.1 : Description of the treatment



Number of session:

1

Total treatment time:

1h00

Description

The ultrasounds of the Beauty Cleaner® deep clean the skin by eliminating impurities and rebalancing pH levels. The vibrations delicately massage the skin; the muscles of the face relax and the skin becomes more elastic.

Effects on the skin

- ▶ Complexion is brighter
- ▶ Impurities are eliminated (epidermis and dermis)
- ▶ Pores are tight and skin surface is smooth
- ▶ Revitalized skin

4.2 : The Protocol



1h00

STEP		TIME	PART OF FACE	TO DO
1		3 MIN	All	Clean the skin (cotton + cream makeup remover + tonic)
2		5 MIN	All	Use the STEAMER (hot steam)
3		10 MIN	1/2 left side (chin, cheeks, around the eyes)	Use the BEAUTY CLEANER® and a mineral water spray ▶ See movements page 27 ▶ The skin must be moistened before each movement of the BEAUTY CLEANER®
		10 MIN	1/2 right side (chin, cheeks, around the eyes)	
		10 MIN	Top of the face (forehead, temples and nose)	
4		10 MIN	All	Use the Ultrasound Probe of the BEAUTY FACE® ▶ Apply the aloe vera gel over the whole face using the Ultrasound Probe
5		5 MIN	All	▶ Remove excess aloe vera gel with cotton and tonic ▶ Perform the facial massage using the day cream
6		3 MIN	All	Manual lymphatic drainage ▶ See movements page 26


4.3 : The Steps



1 3 min Clean the skin 


2 5 min Use the steamer 

3 30 min Beauty Cleaner

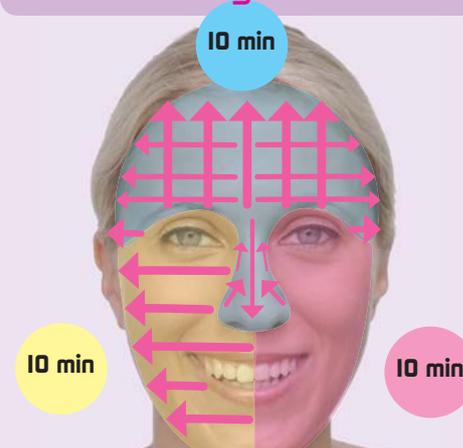

Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.
See details of movement on next page

4 10 min UP Beauty Face

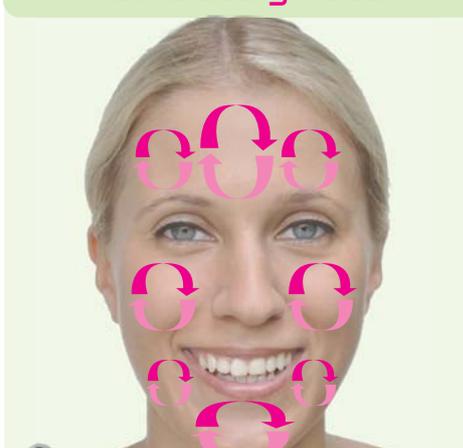


Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

5 5 min Facial Massage



 Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

6 3 min Manual drainage

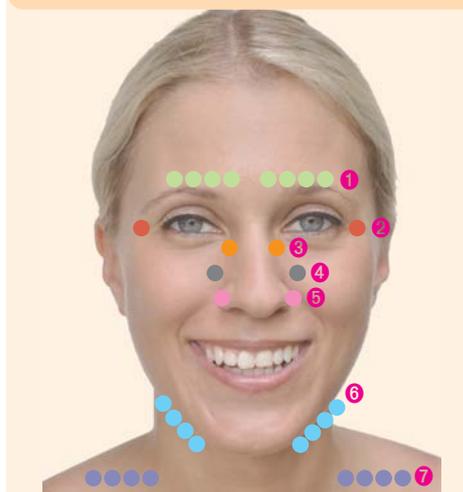

 On each zone (same color) 3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.


4.4 : The Movements

3 30 min Beauty Cleaner®




1

10 min
1/2 left side

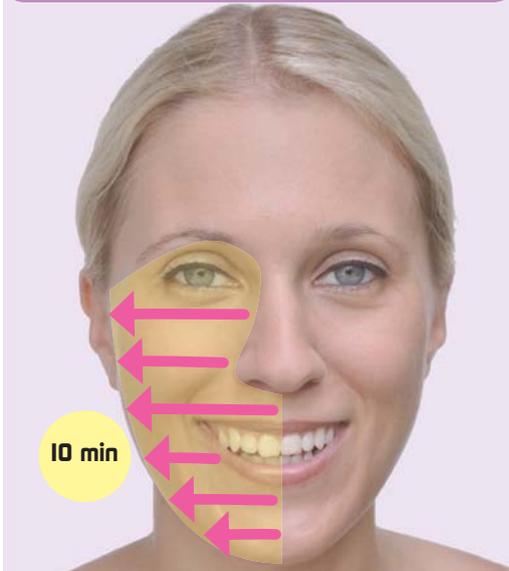


2

10 min
1/2 right side

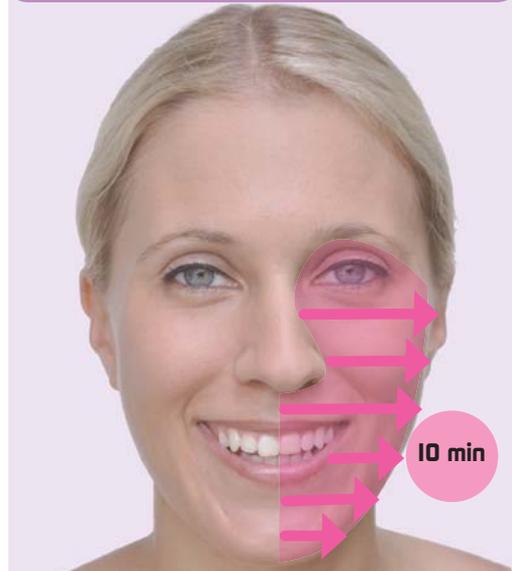

3

10 min
forehead/temple/nose

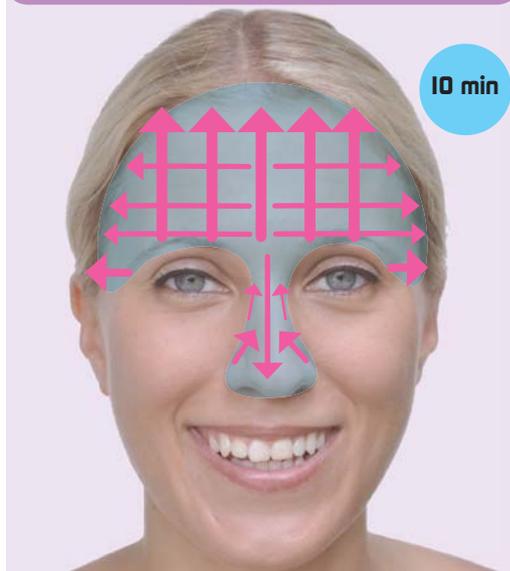


LOW (green):
Sensitive skin and eyes

OR

HIGH (red):
Normal skin

+

PULSE MODE (orange):
For sensitive areas

30 min BEAUTY CLEANER®

- The skin must be moistened before each movement of the Beauty Cleaner.
- For your client's comfort, stop 1 cm before the ear (the vibrations from the ultrasounds can be a nuisance).

Part 5

Treatment Protocols

Debussy® Facial

5.1 : Description of the treatment

5.2 : The Protocol

5.3 : The Steps

5.4 : The Movements



5.1 : Description of the treatment



Number of session:

1

Total treatment time:

1h30

Description

Thanks to the Beauty Cleaner® and to the Beauty Face®, this treatment has a triple effect on the skin: collagen stimulation as well as both deep cleaning and deep hydration of the skin.

Effects on the skin

- Elimination of impurities and toxins (epidermis and dermis)
- Brighter and clearer complexion
- Tightened skin

5.2 :

The Protocol

1 h30



STEP		TIME	PART OF FACE	TO DO
1		3 MIN	All	Clean the skin (cotton + cream makeup remover + tonic)
2		5 MIN	All	Use the STEAMER (hot steam)
3		10 MIN	1/2 left side (from chin to forehead)	Use the BEAUTY CLEANER® and a mineral water spray ▶ The skin must be moistened before each movement of the BEAUTY CLEANER®
		10 MIN	1/2 right side (from chin to forehead)	
4		20 MIN	1/2 left side (from neck to forehead)	Use the BEAUTY FACE® (heating massage) ▶ See movements page 32
		20 MIN	1/2 right side (from neck to forehead)	
5		10 MIN	All	Use the Ultrasound Probe of the BEAUTY FACE® ▶ Apply the aloe vera gel over the whole face using the Ultrasound Probe
6		5 MIN	All	▶ Remove excess aloe vera gel with cotton and tonic ▶ Perform the facial massage using the day cream
7		3 MIN	All	Manual lymphatic drainage ▶ See movements page 31


5.3 : The Steps



1 3 min Clean the skin 


2 5 min Use the steamer 

3 20 min Beauty Cleaner

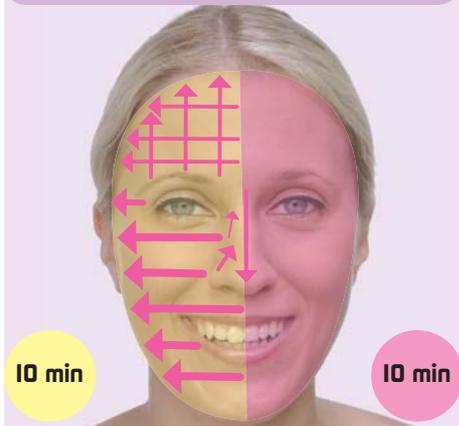


 Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.

4 40 min Facial Tube Beauty Face

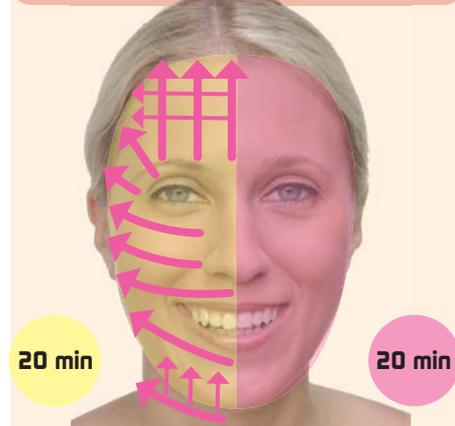

 Normal skin= Hi/Lo: FULL + Power: I
Sensitive skin= FULL + Power: I + Pulse
See details of the movement on next page

5 10 min UP Beauty Face

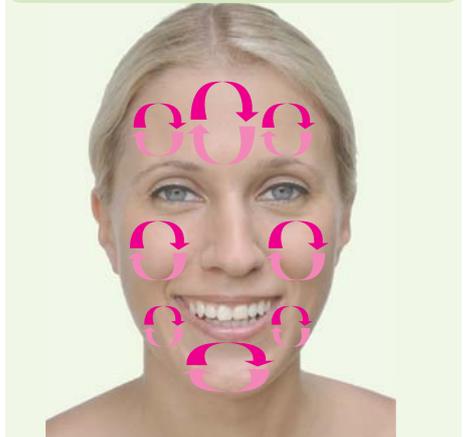


 Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

6 5 min Facial Massage



 Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

7 3 min Manual drainage

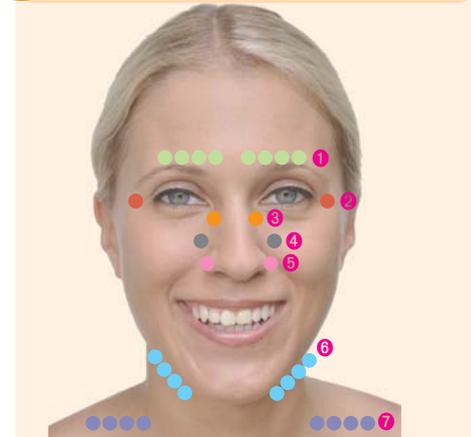

 On each zone (same color) 3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.

40 min
Beauty Face


4




5.4 : The Movements

- ▶ Start from the neck and progress steadily towards the forehead. The heating massage must be performed on the whole surface of the face.
- ▶ One hand uses the tube for the massage while the other hand precedes the movement and stretches the skin in front of the passing tube.
- ▶ Perform long movements starting from the center of the face and stretching towards the ears or the hairline.
- ▶ For a lifting effect, each line rises slightly.
- ▶ When wrinkles are visible, you must stretch the wrinkles with your fingers in front of the passing facial tube.


Neck

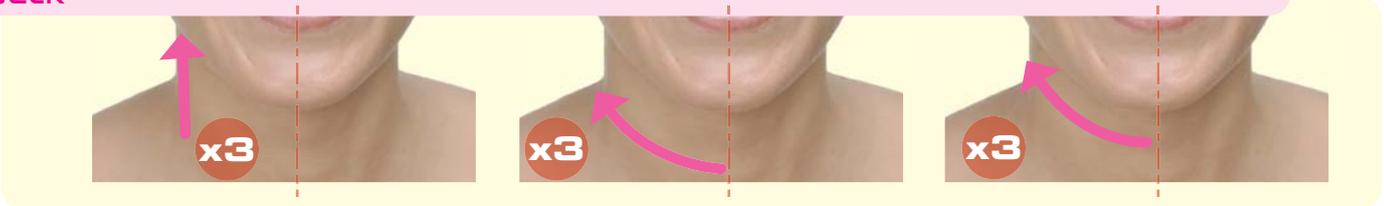

Oval of the face

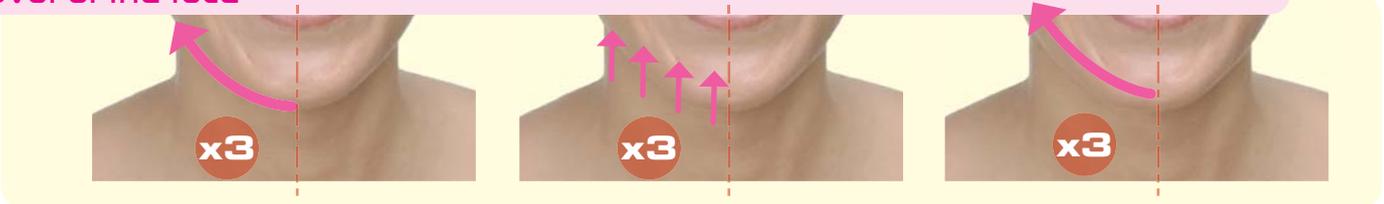

Cheeks

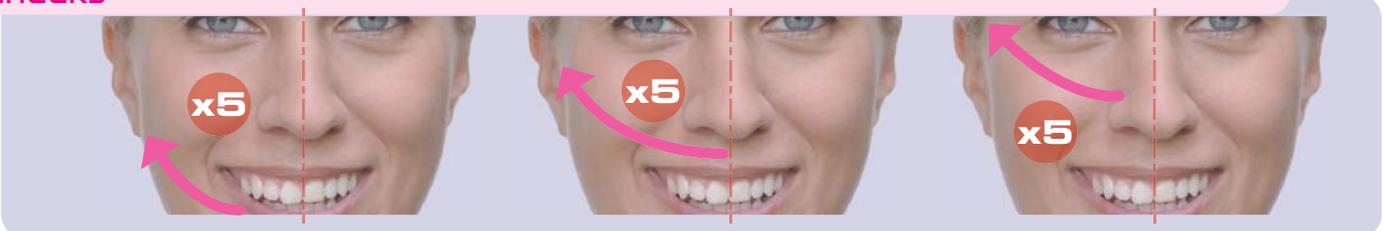

Eyes



Forehead



Part 6

Treatment Protocols

Debussy® Vitality & Debussy® Rejuvenation

6.1 : Description of the treatments

6.2 : The Protocol & Steps
(Session 1)

6.3 : The Protocol & Steps
(Sessions 2,3,4 & 5)

6.4 : The Movements



6.1 : Description of the treatments



	Number of sessions:	Total treatment time:
Debussy® Vitality	3	4h30
Debussy® Rejuvenation	5	7h30

Description

This session treatment provides the skin complete and durable rejuvenation. Your skin will show renewed brightness and greater firmness.

Effects on the skin

- ▶ Oxygenation of tissues and elimination of toxins (epidermis and dermis)
- ▶ Lifting effect, lightening of the skin
- ▶ Bright and clear complexion

Debussy® Vitality

Session 1

1 h30

= Debussy® Facial,

(40 min Beauty Face®)

Session 2

1 h30

= Debussy® Facial,
without Beauty Cleaner®

(1h00 Beauty Face®)

Session 3

1 h30

= Debussy® Facial,
without Beauty Cleaner®

(1h00 Beauty Face®)

Session 4

1 h30

= Debussy® Facial,
without Beauty Cleaner®

(1h00 Beauty Face®)

Session 5

1 h30

= Debussy® Facial,
without Beauty Cleaner®

(1h00 Beauty Face®)

Debussy® Rejuvenation


6.2 : The Protocol and Steps

Session 1
1h30
(same than Debussy® Facial)



STEP		TIME	PART OF FACE	TO DO
1		3 MIN	All	Clean the skin (cotton + cream makeup remover + tonic)
2		5 MIN	All	Use the STEAMER (hot steam)
3		10 MIN	½ left side (from chin to forehead)	Use the BEAUTY CLEANER® and a mineral water spray ➤ The skin must be moistened before each movement of the BEAUTY CLEANER®
		10 MIN	½ right side (from chin to forehead)	
4		20 MIN	½ left side (from neck to forehead)	Use the BEAUTY FACE® (heating massage) ➤ See movements page 37 (same than the Debussy® Facial)
		20 MIN	½ right side (from neck to forehead)	
5		10 MIN	All	Use the Ultrasound Probe of the BEAUTY FACE® ➤ Apply the aloe vera gel over the whole face using the Ultrasound Probe
6		5 MIN	All	➤ Remove excess aloe vera gel with cotton and tonic ➤ Perform the facial massage using the day cream
7		3 MIN	All	Manual lymphatic drainage


3 20 min
Beauty Cleaner

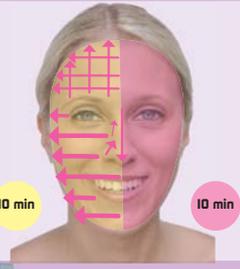

10 min 10 min

Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.


4 40 min
Tube Facial Beauty Face

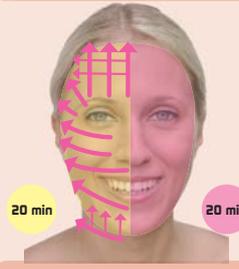

20 min 20 min

Normal skin= Hi/La:
FULL + Power: I
Sensitive skin= FULL + Power :
I + Pulse

5 10 min
UP Beauty Face

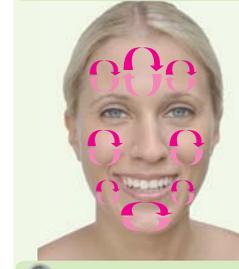


Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

6 5 min
Facial Massage



Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

7 3 min
Manual drainage

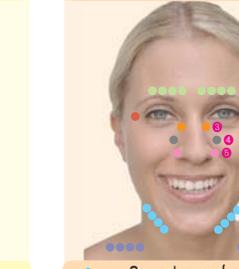

On each zone (same color)
3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.


6.3 : Protocol and Steps



Sessions 2,3,4 & 5
1h30

STEP		TIME	PART OF FACE	TO DO
1		3 MIN	All	Clean the skin (cotton + cream makeup remover + tonic)
2		30 MIN	1/2 left side (from neck to forehead)	Use the BEAUTY FACE® (heating massage) ▶ See movements page 37 (same than the Debussy® Facial)
		30 MIN	1/2 right side (from neck to forehead)	
3		10 MIN	All	Use the Ultrasound Probe of the BEAUTY FACE® ▶ Apply the aloe vera gel over the whole face using the Ultrasound Probe
4		5 MIN	All	▶ Remove excess aloe vera gel with cotton and tonic ▶ Perform the facial massage using the day cream
5		3 MIN	All	Manual lymphatic drainage


2 60 min
Tube Facial Beauty Face

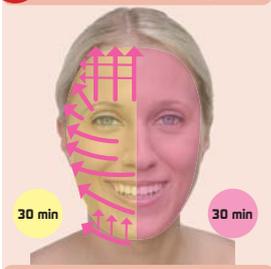

30 min 30 min

 Normal skin= Hi/Lo: FULL + Power: I
Sensitive skin= FULL + Power: I + Pulse

3 10 min
UP Beauty Face

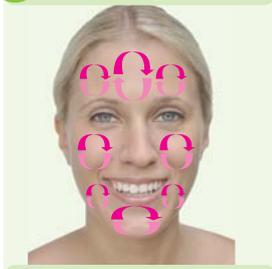


 Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

4 5 min
Facial Massage



 Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

5 3 min
Manual drainage

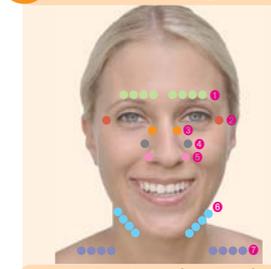


 On each zone (same color) 3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.



6.4 : The Movements

- ▶ Start from the neck and progress steadily towards the forehead. The heating massage must be performed on the whole surface of the face.
- ▶ One hand uses the tube for the massage while the other hand precedes the movement and stretches the skin in front of the passing tube.
- ▶ Perform long movements starting from the center of the face and stretching towards the ears or the hairline.
- ▶ For a lifting effect, each line rises slightly.
- ▶ When wrinkles are visible, you must stretch the wrinkles with your fingers in front of the passing facial tube.

Neck

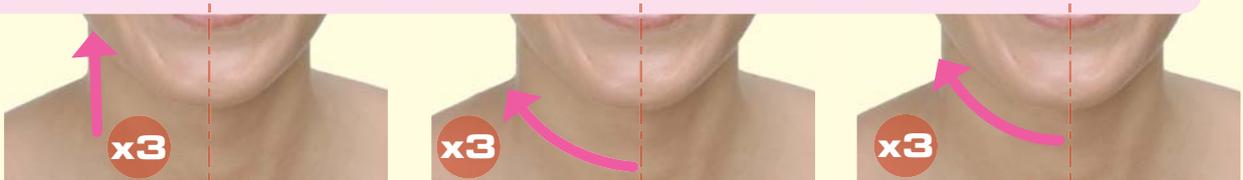

Oval of the face



Cheeks



Eyes



Forehead



Part 7

In the Booth

7.1 : Setting up the booth

7.2 : Keeping your clients informed

7.3 : Your treatment worksheets



7.1 :

Setting up the booth

- ! Your booth must be spotless (clean loor, sheets, pillows...).
- ! The accessories must be cleaned (light, nozzle, probe...).
- ! The facial nozzle for the targeted area is attached to the tube of the BEAUTY FACE®.
- ! The BEAUTY CLEANER® is clean and charged.
- ! The ultrasound probe is plugged into the BEAUTY FACE®.
- ! Creams, aloe vera gel, tonic, makeup remover and cotton are in close proximity to the work space where treatment is being done.
- ! Have at the clients' disposal, eye protection and headband for hair.
- ! Don't forget to put your protection mask on your mouth before start your treatment.

7.2 : Keeping your clients informed

Say...



Present the tools to your client and test heat

«Here, this machine will be used for your treatment. Thanks to its hot massage, infrared and ultrasounds, the collagen will be stimulated, the skin tightened and well nourished. We'll have you test the suction and temperature (test on client's arm). This is what you will feel on your face, it is painless and without danger».

Say...



Step 1, say : «I will start the treatment by cleaning your skin with a cream makeup remover and tonic».



Step 2, say : «then we will use the Steamer to moisten the skin and open your pores. Afterwards I will clean your skin with ultrasounds from the Beauty Cleaner®. During this process, I will spray mineral water over the face».



Step 3, say : «then I will give you a facial massage using infrareds. You will feel some heat, this is normal. If the heat becomes uncomfortable, please tell me»



Step 4, say : «With the help of the Ultrasound Probe, I will help the Aloe Vera penetrate the hypodermis (3cm under dermis) to rehydrate the skin».



Step 5, say : «I'll end the treatment with a facial massage using an Aloe Vera cream and with manual lymphatic drainage».

Debussy® Beauty Cleaner®

1h00

1 3 min
Clean the Skin

2 5 min
Use the Steamer

3 30 min
Beauty Cleaner

Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.

4 10 min
UP Beauty Face

Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

5 5 min
Facial Massage

Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

6 3 min
Manual drainage

On each zone (same color)
3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.

a 10 min
½ left side


b 10 min
½ right side

c 10 min
forehead/temple/nose

1 3 min Clean the skin 

2 5 min Use the steamer 


3 20 min Beauty Cleaner

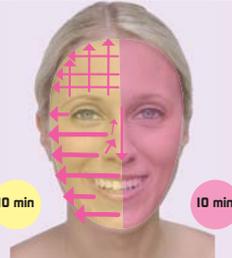

10 min 10 min

Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.

4 40 min Tube Facial Beauty Face

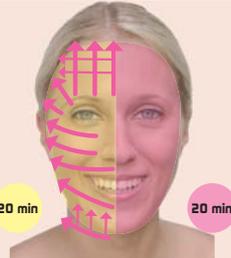

20 min 20 min

Normal skin= Hi/Lo: FULL + Power: I
Sensitive skin= FULL + Power: I + Pulse

5 10 min UP Beauty Face




Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

6 5 min Facial Massage




Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

7 3 min Manual drainage

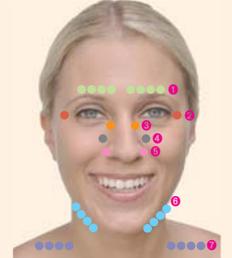

On each zone (same color) 3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.


Neck

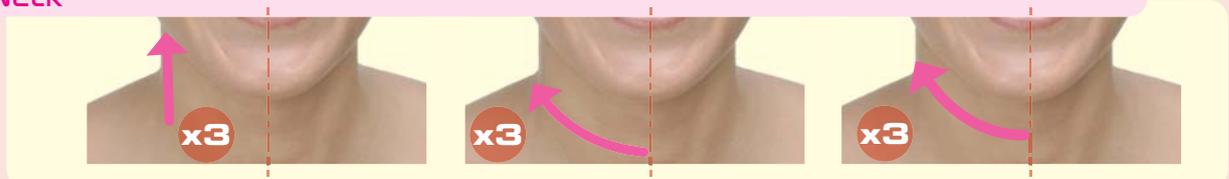

Oval of the face



Cheeks

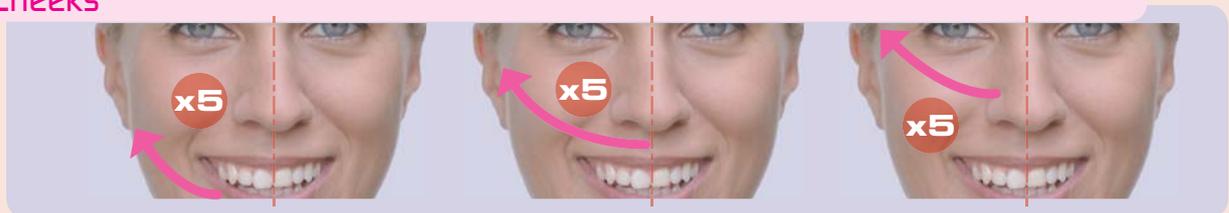

Eyes



Forehead



1 3 min

Clean the skin



2 5 min

Use the steamer



Session 1 : 1h30

3 20 min
Beauty Cleaner

10 min

10 min

Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.

4 40 min (2x20min)
Tube Facial Beauty Face

20 min

20 min

Normal skin= Hi/Lo:
FULL + Power: I
Sensitive skin= FULL + Power: I + Pulse

5 10 min
UP Beauty Face

Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

6 5 min
Facial Massage

Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

7 3 min
Manual drainage

On each zone (same color) 3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.

3 60 min (2x30min)
Tube Facial Beauty Face

Sessions 2 & 3 : 1h30


Neck



Oval of the face



Cheeks

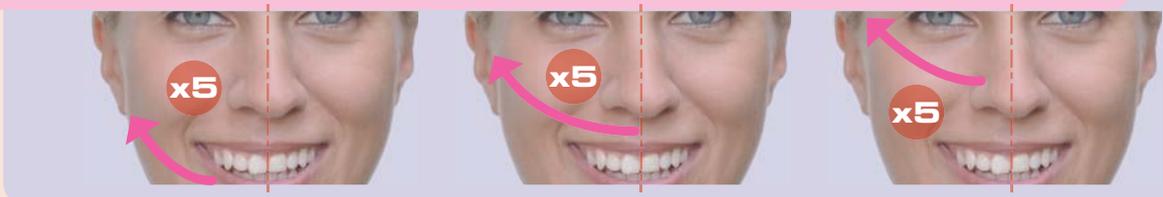

Eyes



Forehead



1 3 min

Clean the skin




2 5 min

Use the steamer




Session 1 : 1h30

3 20 min
Beauty Cleaner

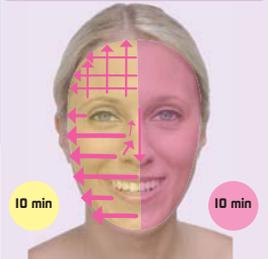


Normal skin= H (High, blue LED)
Moisten the face throughout the treatment.

4 40 min (2x20min)
Tube Facial Beauty Face

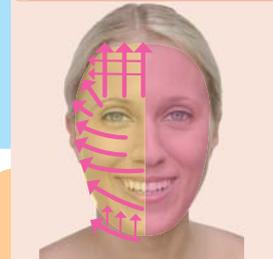

Normal skin= Hi/Lo:
FULL + Power: I
Sensitive skin= FULL + Power: I + Pulse

5 10 min
UP Beauty Face

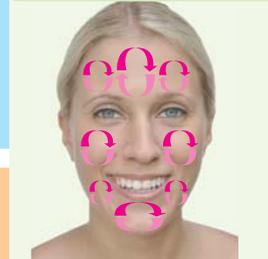


Low to maximum setting
Aloe vera gel must be applied all over the face with the ultrasound probe.

6 5 min
Facial Massage



Remove excess Aloe Vera gel with cotton and tonic.
Apply while massaging the aloe vera day cream.

7 3 min
Manual drainage

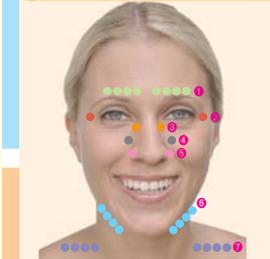

On each zone (same color) 3 medium pressure points upon breathing in.
1 dot = 1 finger, 4 dots = 4 fingers
Use both hands for each zone.

3 60 min (2x30min)
Tube Facial Beauty Face

Sessions 2, 3, 4 & 5 : 1h30

Neck



Oval of the face



Cheeks



Eyes



Forehead



Part 8

Presentation
of the
DEBUSSY®
COSMETIC LINE




Cosmetic Line

Did you know?


For over 5000 years Aloe vera has been used for cosmetic purposes. Aloe vera is 99% water and has a pH level of 4.5 (close to the human skin's). The remaining percentage is composed of over 75 different elements (vitamins, enzymes, minerals, sugars, anthraquinones, lignins, saponins, sterols, salicylic acid and amino acids).

Aloe vera promotes the growth of the skin's fibroblasts cells which produce collagen. It has a firming effect on the skin tissues and slows the appearance of wrinkles.

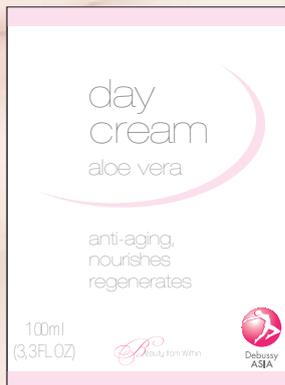


Aloe Vera DAY CREAM

This creamy lotion with ALOE VERA acts in three ways: it nourishes hydrates and slows the skin's aging process. The skin is regenerated, and rediscovers its elasticity, radiance and softness.

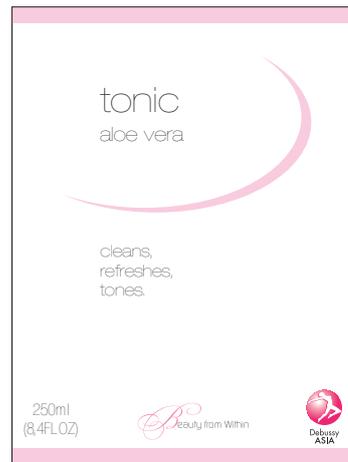


Aloe Vera TONIC

The gentle and nourishing Aloe Vera cleanser hydrates and refreshes the skin. It tones and gives comfort while preserving the skin's natural anti-aging qualities.

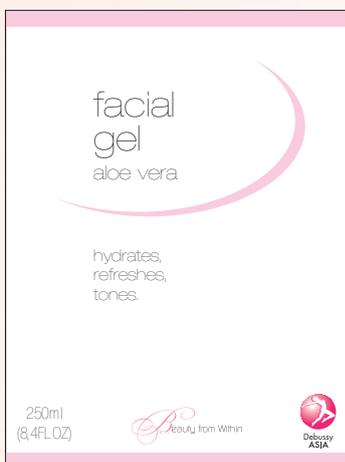

Aloe Vera FACIAL GEL

This refreshing gel reinforces the skin's ability to hydrate itself. It nourishes tones and revitalises. Your skin's tonicity, elasticity and clarity are awakened. It is also advised to use the ALOE VERA Gel after exposure to the sun or directly on sunburns.

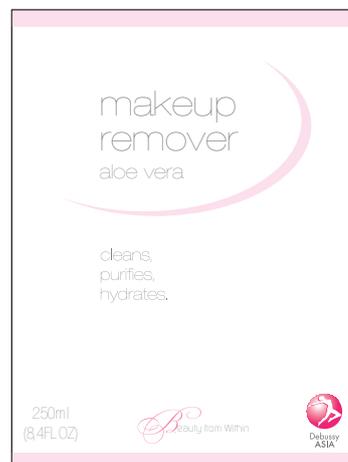

Aloe Vera MAKEUP REMOVER

The makeup remover gently removes makeup and impurities. It respects the skin's natural hydration bringing it more elasticity and freshness.